SmartView Xtend

User Guide

English / Français / Deutsch / Nederlands

This product complies with the requirements of EMC Directive 89/336/EEC.

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European Union Authorised Representative:

Représentant autorisé pour la CEE :
Zugelassene Vertretung in der Europäischen
Gemeinschaft:
Geautoriseerde vertegenwoordiger voor de
Europese Unie:

HumanWare Europe
Russell Smith House
2 Bullmatt Business Centre
Northampton Road, Rushden
Northamptonshire NN10 6AR
United Kingdom
Phone +44 1933 415 800
Fax +44 1933 411 209
Email eu.info@humanware.com

P/N XY947 V1.0

11 Mary Muller Drive, PO Box 3044
Christchurch, New Zealand
phone: +64 3 384 4333 fax: +64 3 384 4933
email: nzinfo@humanware.com
web: www.humanware.com

Declaration of Conformity

(according to ISO/IEC Guide 22 and EN45014)

We: HumanWare
11 Mary Muller Drive
PO Box 3044
Christchurch
New Zealand

Declare that the product:

Product Name: SmartView Xtend

Conforms to the following normative documents:

Safety: BSEN60950-1
BSEN60601-1, IEC601-1, AS/NZS3200.1.0

EMC: BSEN60601-1-2
BSEN55011, CISPR 11
BSEN55022, CISPR 22
BSEN55024, CISPR 24
IEC61000-3-2
IEC61000-3-3
IEC61000-4-2
IEC61000-4-3
IEC61000-4-4
IEC61000-4-5
IEC61000-4-6
IEC61000-4-8
IEC61000-4-11:1994

And complies with the essential requirements of:

Medical Devices Directive	93/42/EEC
Low Voltage Directive	73/23/EEC
EMC Directive	89/336/EEC

Signed:

Russell Watson, Product Manager

Date of Issue:

**Corporate
Headquarters**

HumanWare Limited
11 Mary Muller Drive, P.O. Box 3044
Christchurch, New Zealand
Phone: +64 3 384 4555
Fax: +64 3 384 4933
e-mail: info@humanware.com
Web: www.humanware.com

North America

175 Mason Circle
Concord, CA 94520, USA
Phone: (800) 722 3393
Fax: (925) 681 4630
e-mail: us.sales@humanware.com

1030 René-Lévesque, Drummondville,
Québec, CANADA
Phone: 1-888-723 7273
Fax: (819) 471 4828
e-mail: ca.sales@humanware.com

Europe

Russell Smith House, 2 Bullmatt Business Centre
Northampton Road, Rushden, Northamptonshire NN10 6AR
United Kingdom
Phone: 01933 415800
Fax: 01933 411209
e-mail: eu.sales@humanware.com

Bredaseweg 108a, 4914 NS Oosterhout
The Netherlands
Phone: 31 162 437 910
Fax: 31 162 437 904

Sjouddev. 8, 35246 Växjö
Sweden
Phone: 46 470 764 490
Fax: 46 470 711 088

**Australia
Asia**

Suite 2, 7-11 Railway Street, P O Box 944, Baulkham Hills
NSW 2153, Australia
Phone: (02) 9686 2600
Fax: (02) 9686 2855
e-mail: au.sales@humanware.com

English

SmartView Xtend User Guide

**Copyright © 2006 by HumanWare Limited,
Christchurch, New Zealand. All rights reserved. No part
of this publication may be reproduced or distributed in
any form or by any means, or stored in a database or
retrieval system, without the prior written permission
of HumanWare Limited.**

**Every effort has been made to ensure that the
information contained in this publication is correct at
the time of issue. Neither HumanWare Limited, nor its
agents, assume any responsibility for errors or
omissions. Nor is any liability assumed for damages
resulting from the use of the information contained
herein.**

**11 Mary Muller Drive, P.O. Box 3044
Christchurch, New Zealand
www.humanware.com**

English

Important Safety Instructions

- 1. Follow all warnings, cautions, and instructions marked on the product.**
- 2. The power supply cord is a disconnect device for pluggable equipment. The socket-outlet should be installed near the equipment, and be easily accessible.**
- 3. This product is equipped with a 3-wire grounding type plug as a safety feature. This plug has a third (grounding) pin and can only be plugged into a grounding-type power outlet. If necessary, contact an electrician to replace your power outlet with a grounding-type power outlet.**
- 4. This product should be operated from the type of power indicated on the marking label. If you are not sure of the type of power available, consult your dealer or local power company.**
- 5. Refer all servicing of this product to qualified service personnel. Do not attempt to service it yourself, as opening or removing covers may expose you to dangerous voltage points or other risks.**

Susceptibility to Interference

There may be temporary degradation of the image when the SmartView Xtend is subjected to a strong radio-frequency field, electro-static discharge or

transient electrical noise. An electro-static discharge (caused by static electricity) may cause the screen to go blank. If this happens, switch the unit off and on again.

FCC Warning

This equipment generates and uses radio frequency energy and if not installed and used properly, in accordance with the manufacturer's instruction, may cause interference to radio and television reception. It has been type tested and found to comply with the limits for a Class B digital device in accordance with Subpart B of Part 15 of FCC Rules and Regulations, which are designed to provide reasonable protection against such interference in a residential installation. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:


- 1. Reorient the receiving antenna.**
- 2. Relocate the equipment with respect to the receiver.**
- 3. Move the equipment away from the receiver.**

English

- 4. Plug the equipment into a different outlet so that the equipment and receiver are on different branch circuits.**

If necessary, the user should consult the dealer or an experienced radio/television technician for additional suggestions. The user may find the booklet prepared by the Federal Communication Commission helpful: "How to Identify and Resolve Radio-TV Interference Problems". This booklet is available from the U.S. Government Printing Office, Washington, DC 20402. Stock No. 004-000-00345-4.

Shielded cables and certified Class B peripherals must be used on this product. Using unshielded cables or uncertified peripherals may result in this unit not complying with FCC Rules Part 15. Changes or modifications not expressly approved by the manufacturer could void the user's authority to operate this equipment.

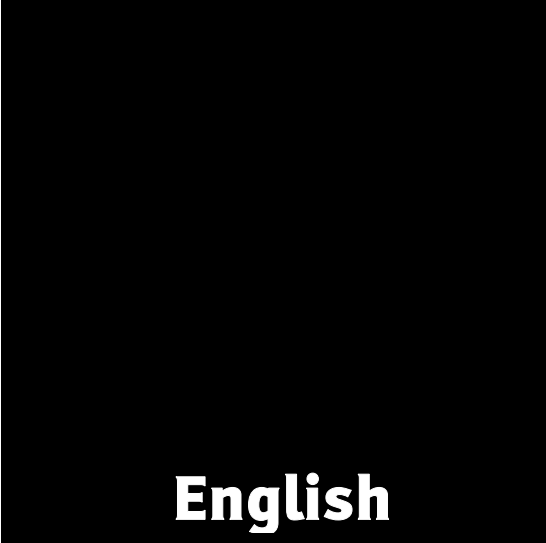English

Contents

Quick Tour - SmartView Xtend 2

Overview of the SmartView Xtend Modules 5

Base Model Controls..... 7

Front Control Panel 7

Margin Stops.....11

Main Menu 12

Xtend Module 1 21

Simple Remote..... 21

Xtend Module 2 29

Keypad 29

Xtend Module 3..... 39

Keypad 39

Preferences Menu 40

Footswitch 41


Maintenance..... 43

Troubleshooting 45


Returning for Service 51

Specifications..... 53

Index 55

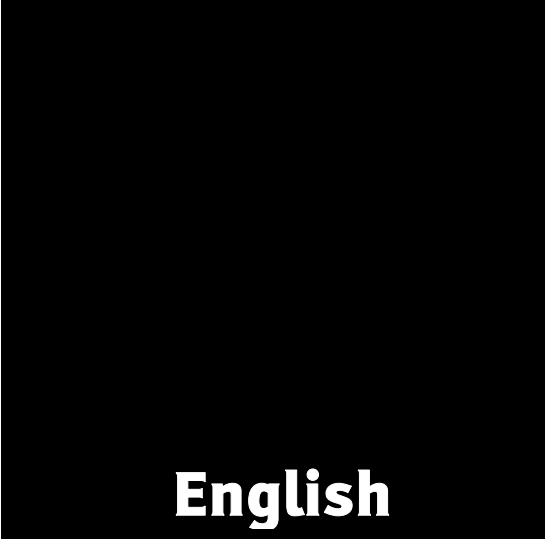
● Quick Tour - SmartView Xtend

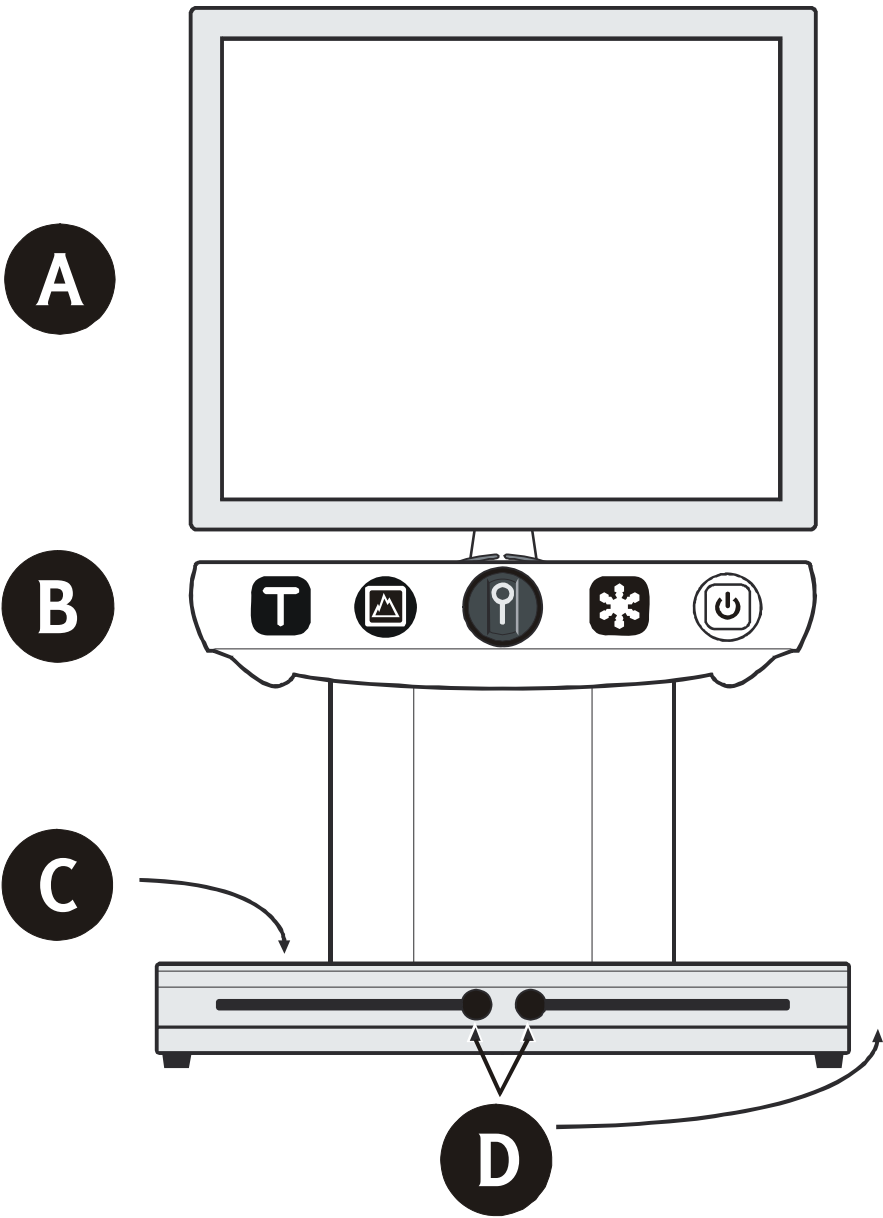
- | | | | |
|----------|---------------------|----------|---------------|
| A | Screen | C | Reading Table |
| B | Front Control Panel | D | Margin Stops |

English

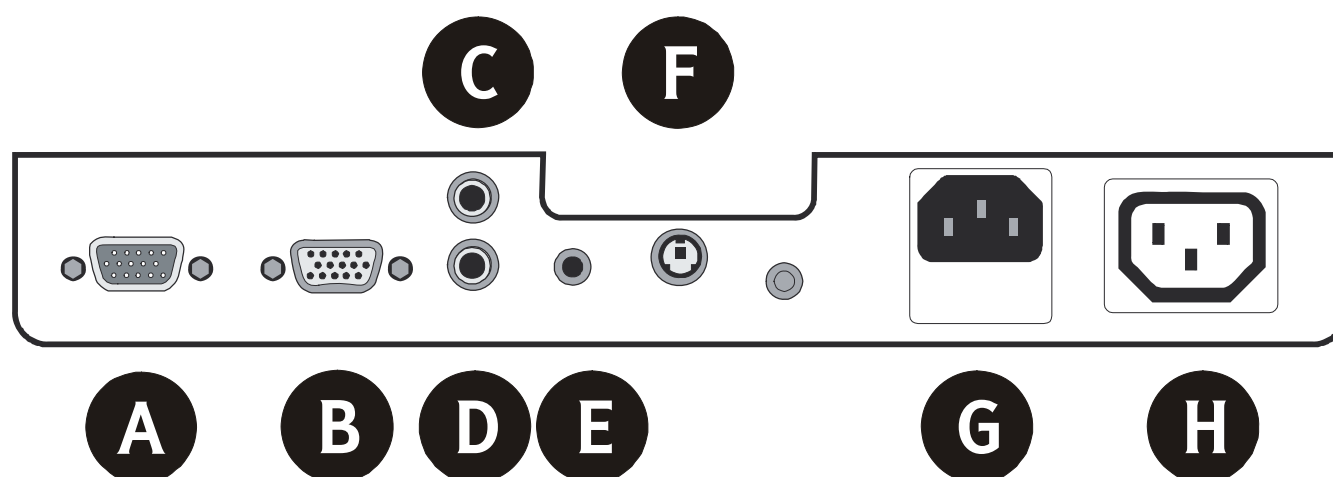- A** PC Screen
- B** SmartView Screen
- C** Video Camera Input/TV In
- D** TV output (if used instead of a VGA Monitor)
- E** Footswitch
- F** Simple Remote or Keypad
- G** Power Outlet (e.g., to Monitor)
- H** Power Input

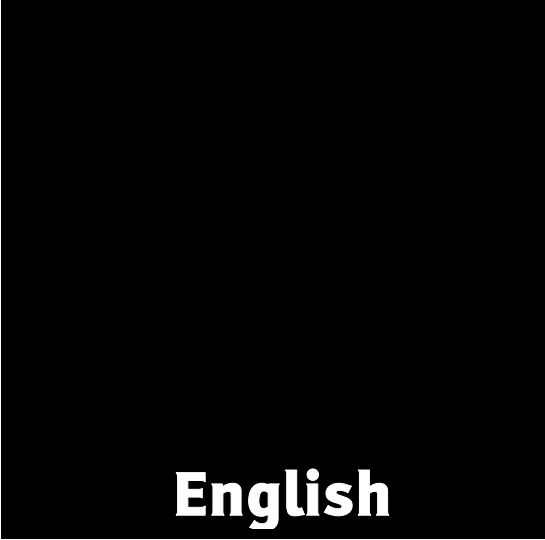● **Overview of the SmartView Xtend Modules**

► **Module 1**

In addition to the Front Control Panel, Module 1 includes a Simple Remote with 3 extra functions:

- A Page Locator to highlight your place on a page.
- A Preset function to program your preferred zoom setting into the machine.
- A Lines and Blinds function to set reading guides.

► **Module 2**

Module 2 includes a Keypad with the same functions as the Simple Remote, along with 3 extra functions:

- An onscreen clock and calendar.
- A calculator.
- A Preferences menu that allows you to change more settings.

► **Module 3**


Module 3 includes the Keypad above with all the same functions, as well as an option to connect a Personal Computer (PC) or an external video camera to the SmartView Xtend.

English

You can upgrade your SmartView to another module at any time. Contact your local distributor for details.

● Base Model Controls

■ Front Control Panel

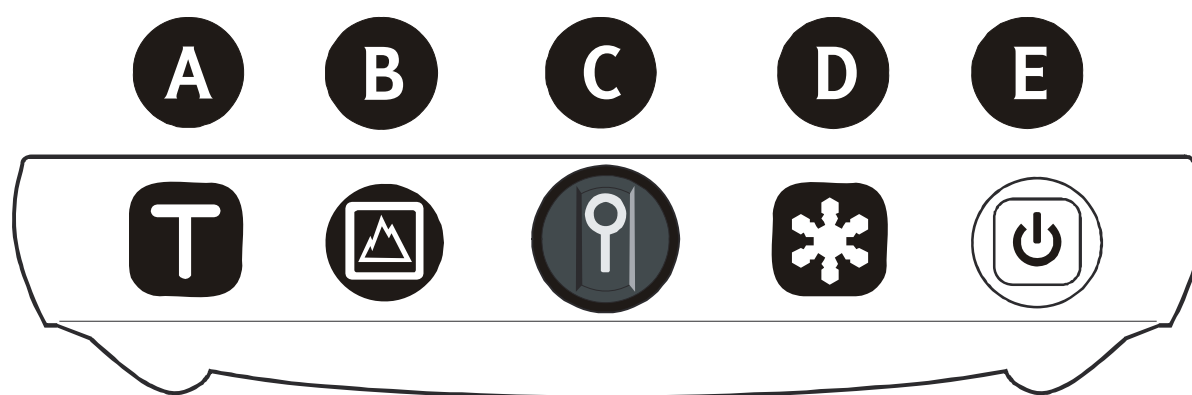- A** Text Button
- B** Picture Button
- C** Zoom Knob
- D** Brightness Button
- E** Power Button

English

► **Power Button**

- To turn your SmartView Xtend on and off, press the Power button.

You will know that the machine is on when the lights over the reading table illuminate and the green light in the Power button glows.

► **Text Button**

Use this button if you want to read text.

False Colour Combinations

There are several false colour combinations available in Text mode.

- To view the 18 combinations of foreground and background colours, press the Text button repeatedly.

The display cycles through the first 9 combinations in the following order:

- Black text on a white background.
- Blue text on a white background.
- Blue text on a yellow background.
- Black text on a yellow background.
- Black text on a purple background.
- Black text on a light blue background.
- Black text on an orange background.

- Black text on a green background.
- Black text on a grey background.

The display then cycles through the opposites of the combinations above.

① Handy Hint

If you prefer to look at text on a white background, but find the glare too much with a particular document, try using black text on a grey background.

Preferred False Colour Combination

If you use a particular colour combination often, you can change the display so that the SmartView Xtend only cycles through:

- Your preferred colour combination.
- White on black.
- Black on white.
- To do this:
 1. Select your preferred colour combination.
 2. Hold down the Text button until you hear a beep.

Note: SmartView Xtend retains all of your settings, even when it is turned off and then on again.

- To restore the full set of 18 colour options, hold down the Text button and the Picture button together until you hear a beep.

English

► **Picture Button**

Use this button if you want to view photos or pictures.

► **Zoom Knob**

Use this knob to adjust the size of your displayed image.

- To enlarge the image, turn the Zoom knob clockwise.
- To reduce the image, turn the Zoom knob anti-clockwise.

► **Brightness Button**

Use this button to increase or decrease the brightness of the display.

- To adjust the brightness:
 1. Press the Brightness button. A slide bar appears on the display.
 2. Move the slide bar between 0%-100% by turning the Zoom knob clockwise to increase the brightness, and anti-clockwise to decrease it.
 3. Press the Text or Picture button to return to the normal display.

Note: You can also adjust the brightness levels within the menu. See Menu Controls/Menu Items/Brightness (on page 14) for more details.

① Handy Hints

Reducing Glare

If you prefer to look at text on a white background but find the glare too much with a particular document, reduce the brightness. You can increase it again later if necessary.

Setting Brightness for Text or Photos

When you change the brightness of the display, the adjustment you make only occurs within the current mode. If you change the brightness level in Text Mode, the level in Picture Mode remains unchanged.

This can be useful for several reasons. You may find that a high level of brightness causes too much glare in Text Mode, for example, but that the same level works well when you are viewing pictures.

■ Margin Stops

Use the Margin Stops to reduce the movement range of the reading table for a particular document. This can prevent your reading material from disappearing off the edge of the screen as you read.

- To set margins:
 1. Line up the left-hand margin of your reading material with the screen.

English

2. Unscrew the left-hand Margin Stop.
3. Slide the left-hand Margin Stop as far to the right as possible and tighten it.
4. Line up the right-hand margin of the reading material with the screen.
5. Unscrew the right-hand Margin Stop.
6. Slide the Margin Stop as far to the left as possible and tighten it.

■ Main Menu

The Main Menu allows you to adjust a number of different settings.

- To enter the Main Menu, hold down the Brightness button until the menu appears.

The Main Menu includes the following items:

- Brightness.
- Contrast.
- Colour.
- Normal/Enhanced Colour.
- Output Display Frequency.
- Beeps.
- Language.
- Region.
- System Information.
- Preferences (Modules 2 and 3 only).

Note: If you do not touch any buttons or knobs for 20 seconds when a menu is displayed, the display reverts to normal operation.

► **Menu Controls**

When you are using the menu, the functions of the Front Control Panel are as follows:

Brightness Button

- Use the Brightness button to:
 - Display the Main Menu.
 - Select a menu item.

Zoom Knob

- Use the Zoom knob to:
 - Move up and down within the menu. (Turn the knob clockwise to move down the menu items and anti-clockwise to move up.)
 - Slide between 0% and 100% when a slide bar is displayed - for example, when adjusting the brightness or contrast.

Picture Button

- Use the Picture button to:
 - Move up through the levels of the menu system.

English

For example, after adjusting the contrast, you can return to the Main Menu by pressing the Picture button.

- Exit the Main Menu.

Text Button

- Use the Text button to exit the menu at any time.

► Menu Items

Brightness

Use this menu item to increase or decrease the brightness of the display. You can do this independently in Text and Picture modes.

- To adjust the brightness via the menu:
 1. Navigate to Brightness by turning the Zoom knob.
 2. Select Brightness by pressing the Brightness button.
 3. Move the Brightness slide bar between 0%-100% by turning the Zoom knob.
- To return to the Main Menu, press the Picture button.
- To exit the menu and return to the normal display, press the Text button.

Contrast

Use this menu item to adjust the contrast of the display. You can do this independently in both Picture or Text mode.

- To adjust the contrast:
 1. Navigate to Contrast by turning the Zoom knob.
 2. Select Contrast by pressing the Brightness button.
 3. Move the Contrast slide bar between 0%-100% by turning the Zoom knob.
- To return to the Main Menu, press the Picture button.
- To exit the menu and return to the normal display, press the Text button.

Colour

Use this menu item to adjust the colour of a picture or photo from fully enhanced colour to grey-scale.

- To adjust the colour of an image:
 1. Navigate to Colour by turning the Zoom knob.
 2. Select Colour by pressing the Brightness button.
 3. Move the Colour slide bar between 0%-100% by turning the Zoom knob.

English

- To return to the Main Menu, press the Picture button.
- To exit the menu and return to the normal display, press the Text button.

Normal/Enhanced Colour


Use this menu item to view a picture in either realistic colour or moderately enhanced colour.

- To enhance the colour of an image:
 1. Navigate to Normal/Enhanced colour by turning the Zoom knob.
 2. Switch between Normal and Enhanced colour by pressing the Brightness button.
- To exit the menu, press the Text or Picture button.

Output Display Frequency

Use this menu item to change the VGA display frequency from 60 Hz to 75 Hz.

- To adjust the display frequency:
 1. Navigate to Output Display Frequency by turning the Zoom knob.
 2. Switch between 60 Hz and 75 Hz by pressing the Brightness button.
- To exit the menu, press the Text or Picture button.

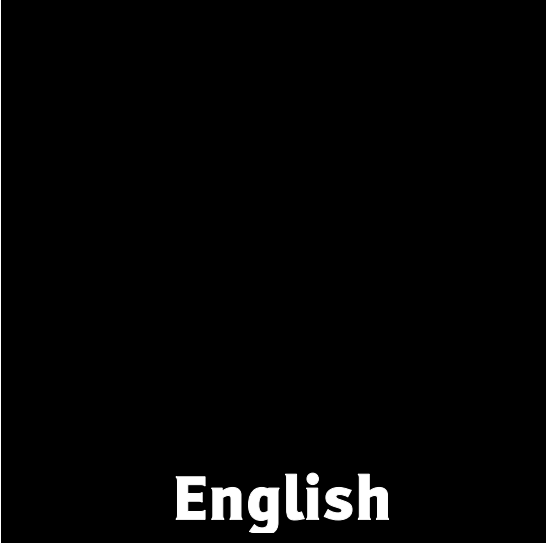English

Language

Use this menu item to set the language of the onscreen display.

- To set the language:
 1. Navigate to Language by turning the Zoom knob.
 2. View the supported languages by pressing the Brightness button.

The display cycles through English, Español, Nederlands, Italiano, Français, Deutsch, Svenska, Português, Norsk, English US.

To select a language, keep pressing the Brightness button until you reach the language that you want. All the menu items then change accordingly.

- To exit the menu, press the Text or Picture button.

Region

Use this menu item to set the region, which determines the contents of the System Information screen. See Menu Items/System Information (on page 18).

- To set the region:
 1. Navigate to Region by turning the Zoom knob.
 2. View the supported regions by pressing the Brightness button.

English

The display cycles through United States, Canada, United Kingdom, New Zealand, Australia & Asia, Netherlands, Sweden, and none.

To select a region, keep pressing the Brightness button until you reach the region that you want. The support contact details on the System Information screen automatically change to match the selected region.

- To exit the menu, press the Text or Picture button.

Beeps

Use this menu item to switch the audible beeps on or off.

- To adjust the Beeps setting:
 1. Navigate to Beeps by turning the Zoom knob.
 2. Switch between Beeps On and Beeps Off by pressing the Brightness button.
- To exit the menu, press the Text or Picture button.

System Information

Use this menu item to view the system configuration information.

- To view the System Information screen:
 1. Navigate to System Information by turning the Zoom knob.

2. Select System Information by pressing the Brightness button.

- **To return to the Main Menu, press the Picture button.**
- **To exit the menu and return to the normal display, press the Text button.**

Preferences

This menu item is available with upgrade Modules 2 and 3.

See Xtend Module 2/Preferences Menu (on page 35).

● Xtend Module 1

Module 1 provides all the functionality of the Front Control Panel as well as a Simple Remote with several extra functions.

■ Simple Remote

Units with upgrade Module 1 include the Simple Remote.

Note: You cannot use the Simple Remote and the Keypad at the same time.


The following Simple Remote buttons work in the same way as the corresponding buttons on the Front Control Panel. Refer to:

- Text Button (on page 8).
- Picture Button (on page 10).
- Zoom Knob (on page 10).
- Brightness Button (on page 10).

The remaining buttons do not appear on the Front Control Panel and have their own unique functions.

English

Xtend Module 1

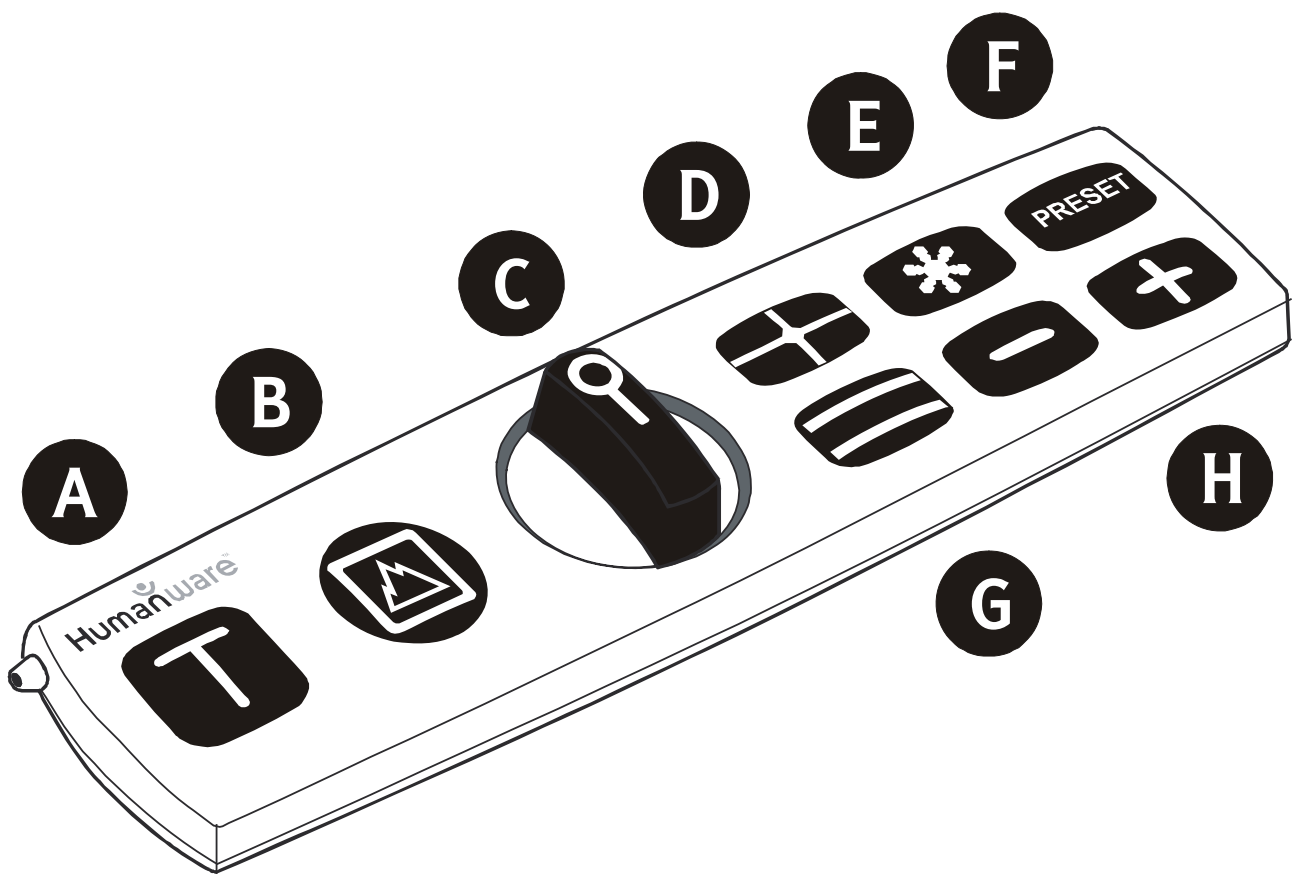- A** Text Button
- B** Picture Button
- C** Zoom Knob
- D** Page Locator
- E** Brightness Button
- F** Preset Button

G Lines/Blind Button

H Plus & Minus Buttons

► Page Locator Button

Use this button to find your position on a page.

- To find out your location on your reading material, press the Page Locator button.

A red dot shines down onto the reading material, and a red light is displayed at the corresponding point in the centre of the screen.

- To turn off the red dot, press the Page Locator button again.

① Handy Hint

If you want to make sure that you sign your name in the correct place when signing a document, you can use the Page Locator as follows:

1. Press the Page Locator button. A red dot appears on the reading table.
2. While looking at the screen, position the document so that the red light on the screen is aligned with where you need to start writing.

English

3. Now look at the reading table and position your pen on the red dot. Look at the screen again. Your pen should be in the correct place.

4. Begin writing.

► **Preset Button**

Use this button to preset a zoom setting of your choice. For example, if there is a particular level of magnification that you like to use often, you can store the level in the machine.

- To store a zoom setting, hold down the Preset button until you hear a beep.
- To select the stored zoom setting when you are in normal zoom, press and release the Preset button.

① **Handy Hint**

Preset the magnification for documents that you read regularly (such as the newspaper).

Then whenever you change the magnification to read other documents, you can simply push the Preset button to restore the magnification to read the newspaper again.

► **Lines and Blinds Button**

Use this button to set reading guides.

- To view the reading guide options, press the Lines and Blinds button repeatedly.

The display cycles through:

- Horizontal lines to assist with general reading.
 - Vertical lines to assist with reading a column of text.
 - Horizontal blinds that blank out the top and bottom of the display to remove adjacent lines of text.
 - Vertical blinds that blank out the sides of the display to assist with reading a column.
 - No lines or blinds.
- To change the location or size of the lines and blinds, press the + and - buttons. See below.

► **Plus (+) and Minus (-) Buttons**

Use these buttons to set the spacing of lines and blinds, or to shift them across the screen, as described below:

- To widen the spacing of lines or blinds, press the + button repeatedly until you achieve the desired amount of separation.
- To narrow the spacing of lines or blinds, press the - button repeatedly until the lines or blinds are as close as you want them.

- To move a pair of lines or blinds up or down, hold down both the + and - buttons until a beep sounds. Then press + to move up, or - to move down.
To restore normal operation, hold down both the + and - buttons again for several seconds.

① Handy Hints

Viewing cluttered pictures

If you are viewing a picture with a large amount of activity or clutter in the background, you can use the blinds to conceal the background, as follows:

1. Press the Lines and Blinds button until you reach the appropriate blinds option.
2. Adjust the blinds so that they are similar in size to the foreground area by pressing + and -. See Plus (+) and Minus (-) Buttons (on page 25).
3. Move the blinds over the foreground by pressing + and -.

The foreground items should be easier to see now that the background is obscured.

Reading

If you are finding it difficult to navigate with the reading table, you can use the horizontal blinds to guide your reading, as follows:

- Press the Lines and Blinds button until you reach the horizontal blinds option, then use the + and - buttons to:
 1. Set the blinds to a width of 1-3 lines of text.
 2. Move the blinds to the top of the document.
 3. Move the blinds downwards as you read, instead of moving the reading table.

See Plus (+) and Minus (-) Buttons (on page 25).

● Xtend Module 2

Module 2 provides all the functionality of Module 1 as well as a Keypad with several extra functions.

■ Keypad

Module 2 and Module 3 include the Keypad.


Note: You cannot use the Keypad and the Simple Remote at the same time.

The Keypad provides access to:

- An onscreen clock and calendar.
- A calculator.
- Reading guides (lines and blinds).
- The onscreen menu.
- A Preset function (to adjust your zoom).
- A Page Locator.

English

A B C D E

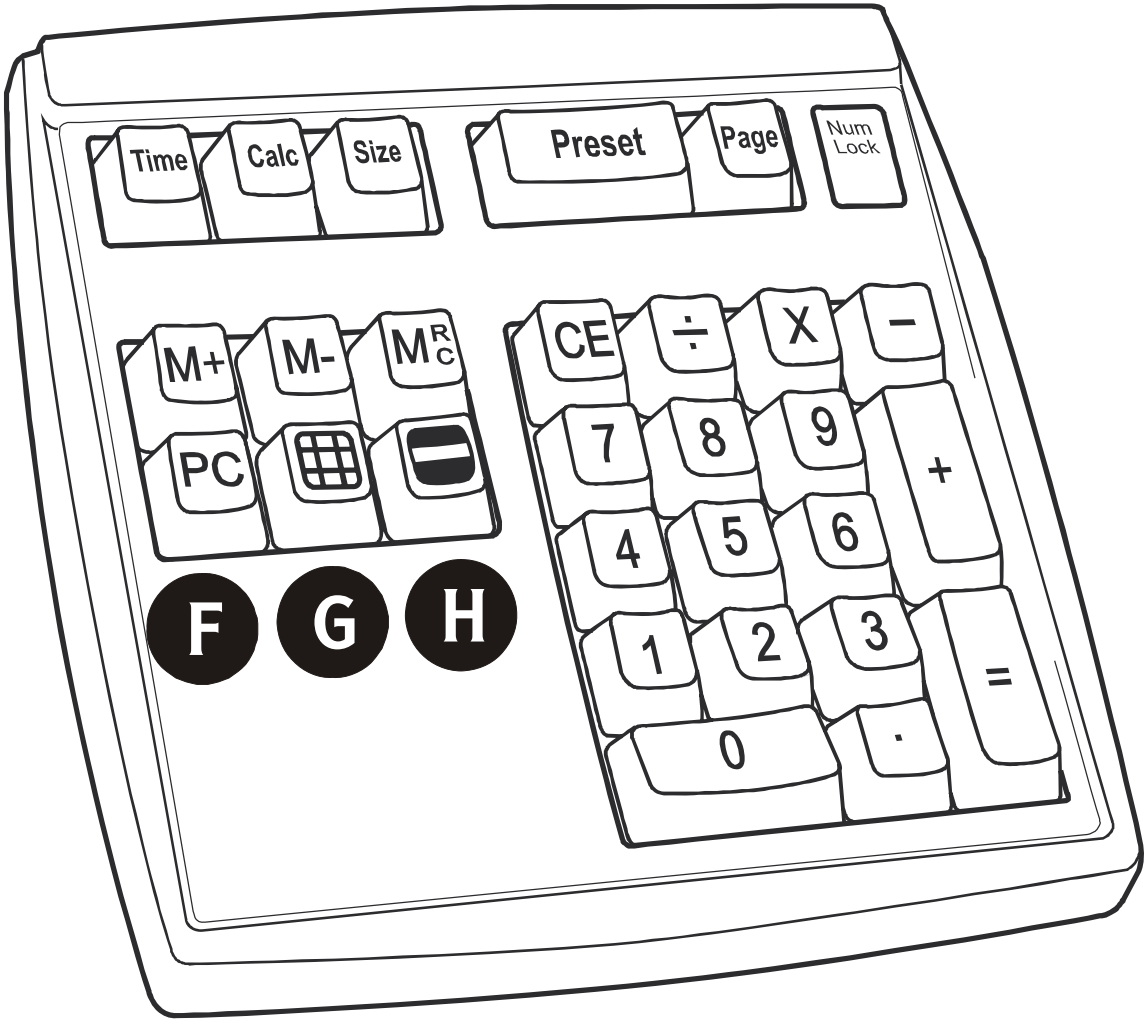Xtend Module 2

- A** Time Key
- B** Calc Key
- C** Size Key
- D** Preset Key
- E** Page Key
- F** PC Key
- G** Lines Key
- H** Blinds Key

► **Time Key**

Use this key to display the time and date in the corner of the screen.

- To turn the time and date display on or off, press the Time key.

Setting the Time and Date

1. To highlight the hour digits, press the Time key for at least one second.
2. To adjust the hour digits, press the ↑ and ↓ keys.
3. To move through the hours, minutes, seconds, and date settings from left to right, press the → key.
4. To adjust each setting (i.e. minutes, seconds, and date), press the ↑ or ↓ keys.
5. To move back through the settings from right to left, press the ← key.
For example, if you accidentally move from Hours to Minutes before you have adjusted the Hour values, you can press ← to go back.
6. To set the time and date and exit the setting procedure, press =.

► Calc Key

Use this key to display the calculator in the top left corner of the screen.

- To turn the calculator on or off, press the Calc key.

Using the Calculator

- To operate the calculator, use the function keys on the Keypad (+, -, *, and /). You can enter up to nine digits.
- To clear the display, press C.
- To add the displayed number to the memory, press M+.
- To subtract the displayed number from the memory, press M-.
- To display the contents of the memory, press MRC (Memory Recall/Clear) once.
- To clear the memory, press MRC twice.

Note: The Time, Date, and Calculator are always displayed respectively in vertical order. They also always appear in the opposite of the current colour combination, although this is only obvious in Text Mode.

► Size Key

Use this key to adjust the size of the Time, Date, and Calculator displays.

- To view the available size options, press the Size key repeatedly. The display cycles through five different sizes of onscreen text.

► **Preset Key**

This key functions in the same way as the Preset button on the Simple Remote. See Simple Remote/Preset Button (on page 24).

► **Page Key**

This key functions in the same way as the Page Locator button on the Simple Remote. See Simple Remote/Page Locator Button (on page 23).

► **PC Key**

This key is only usable in Module 3.

► **Lines Key**

Use this key to set line markers as reading guides.

- To view the available options, press the Lines key repeatedly.

The display cycles through:

- Horizontal lines to assist with general reading.
- Vertical lines to assist with reading a column of text.
- No lines.

► **Blinds Key**

Use this key to set blinds as reading guides.

- To view the available options, press the Blinds key repeatedly.

The display cycles through:

- Horizontal blinds that blank out the top and bottom of the display to remove adjacent lines of text.
- Vertical blinds that blank out the sides of the display to assist with reading a column.
- No blinds.

► **Arrow Keys**

Use the Arrow Keys to:

- Move lines and blinds up or down, using ↑ and ↓.
- Increase or decrease the separation of lines and blinds, using ← and →.
- Change the position of the Time, Date, and Calculator displays on the screen.
- Move up and down within the onscreen menu. Make sure that the full-screen camera image is displayed and the lines or blinds are turned off.

► **Preferences Menu**

The Preferences menu is available with upgrade Modules 2 and 3.

To view the Preferences menu:

1. Open the Main Menu by holding down the Brightness button.
2. Navigate to Preferences by turning the Zoom knob.
3. Select Preferences by pressing the Brightness button.
4. Move through the items in the Preferences menu by turning the Zoom knob on the Front Control Panel or pressing the ↑ and ↓ keys on the Keypad.
5. Cycle through the options for each item by pressing the Brightness button on the Front Control Panel or the ← or → keys on the Keypad.

The Preferences menu includes:

- A separator toggle for the calculator so that you can choose either a fullstop or comma for the calculator's decimal point.
- An operator toggle for the calculator, to switch the function key display (+, -, *, /, +) On and Off.
- A clock toggle, to switch between the 12 and 24 hour clock.
- A time out function, to set the time delay before the lamps turn off when a PC is attached.

Whenever you display an image from a computer or an external video camera, the lamps turn off automatically if you have not displayed the computer or camera image for a certain period of time. You can set this period to be anywhere

between 2 seconds and 10 minutes, or you can set the machine to the Lamps Always On mode.

- A Mode text toggle that you can switch On and Off.

If this option is turned on, the machine will display a short message when switching between Text and Photo modes.

Note: The machine remembers all settings, even when you have switched it off.

● Xtend Module 3

Module 3 provides all the functionality of Module 2 as well as the ability to display images from a Personal Computer (PC) and external video camera.

■ Keypad

Units with Module 2 and Module 3 include the Keypad.
See Xtend Module 2/Keypad (on page 29).

► PC Key

With upgrade Module 3, the PC key on the Keypad allows you to display the contents of a PC screen or external video camera on your SmartView Xtend.

Note: You can only view external video camera images when using the Footswitch. See Xtend Module 3/Footswitch (on page 41).

The Display Options

- To view the display options, connect a PC to your SmartView Xtend using the cable provided, then hold down the PC key.

The display cycles through:

- The SmartView Xtend screen.

English

- A split screen, with one side SmartView Xtend and one side PC.
- Full PC screen.
- To adjust the split screen:
 - Press the Lines key on the Keypad to switch between a horizontal and vertical split.
 - Press the Blinds key on the Keypad to switch between top and bottom and left and right.
 - Use the appropriate Arrow Keys on the Keypad to move the split left, right, up, or down.
 - Press C on the Keypad at any time for a 50% split in vertical and horizontal modes. This is the factory setting.

Note: Split-screen mode is only available when a VGA monitor is attached to the VGA output. If there is no VGA monitor attached, you will hear a warning beep when you press the PC key.

■ Preferences Menu

The Preferences menu is available with upgrade Modules 2 and 3. See Xtend Module 2/Preferences Menu (on page 35).

■ Footswitch

The Footswitch is included with Module 3. You can use it to display the contents of a PC screen and/or external camera on your SmartView Xtend.

Display Options

- To view the display options, connect a PC and/or external camera to your SmartView Xtend using the cable provided, then press the Footswitch repeatedly.

If both a PC and external camera are connected, the display cycles through:

- SmartView Xtend screen.
- Split screen, with one side SmartView Xtend and one side PC.
- Full PC screen.
- Split screen, with one side SmartView Xtend and one side external video camera.
- Full external video camera screen.
- To adjust the split screen:
 - Press the Lines key on the Keypad to switch between a horizontal and vertical split.
 - Press the Blinds key on the Keypad to switch between top and bottom and left and right.
 - Use the appropriate Arrow Keys on the Keypad to move the split left, right, up, or down.

English

- Press C on the Keypad at any time for a 50% split in vertical and horizontal modes. This is the factory setting.

Note: Split-screen mode is only available when a VGA monitor is attached to the VGA output. If there is no VGA monitor attached, you will hear a warning when you press the Footswitch.

English

● **Maintenance**

Wipe your SmartView Xtend with a clean, soft cloth.

● **Troubleshooting**

▶ **If the Keypad is not responding:**

- Check that the Keypad plug is properly inserted into the correct socket on the rear panel.

▶ **If the Simple Remote is not responding:**

- Check that the Simple Remote plug is properly inserted into the correct socket on the rear panel.

▶ **If the Footswitch is not responding:**

- Check that the Footswitch plug is properly inserted into the correct socket on the rear panel.

▶ **If the screen is blank:**

- Check that the AC cables are plugged in properly.
- Check that the power is switched on at the wall socket.
- Check that the green power light is glowing.
- Check that all cables are properly plugged in.
- Check that the monitor is switched on.
- Check that the Brightness and Contrast levels are turned up.
- Turn the Zoom knob anti-clockwise for minimum magnification.

► **If there is no computer display when a PC is connected:**

- Make sure that all connectors are properly plugged in. Particularly make sure that the computer cable is connected to the Xtend unit and not directly to the monitor.
- Make sure that your machine has been upgraded to Module 3.
- Press the PC Key.
- Make sure that your computer is on and working by connecting the PC directly to the monitor.
- Seek help from a qualified professional if all of the above suggestions fail.

► **If the image is out of focus:**

- Try moving your reading material off the reading table and then putting it back again.
- Clean the monitor screen with a glass cleaning solution.

► **If the text looks grey or low in contrast:**

- Press the Text button to set the unit in High Contrast mode.
- Check the Brightness, Contrast, and Colour settings on the menu. Colour should be set to around 100%, and Brightness and Contrast should be set to at least 50%.
- Check the Brightness and Contrast on your monitor. (**Caution required:** seek help from a qualified person.)
- Check that the lights over the reading table are working.

► **If the colours are difficult to detect:**

- Check that the colour setting in the menu is 100%.
- Check that the lights over the reading table are working.

► **If the contrast is low:**

- Check the Contrast setting in the menu. Remember that there are different settings for Picture and Text modes.

English

► **If the brightness is low:**

- Check the Brightness setting in the menu.
Remember that there are different settings for Picture and Text modes.
- Check that both lights over the reading table are working.

► **If there is a lot of glare:**

- Reduce the Brightness setting in the menu.
- Reduce the external lighting effects.

► **If the reading table does not move over the expected range:**

- Check the Margin Stops.
- Check for obstructions under the reading table.

► **If the Page Locator light does not appear:**

- Remove any paper from the reading table and then see if the Page Locator light is visible against the surface of the table.

Certain colours on your reading material may camouflage the colour of the Page Locator light.

► **If there is a limited number of false colours:**

- Hold down the Text and Picture buttons together for five seconds until you hear a beep.

► **If you are stuck in the menu:**

- Press the Text button.
- Press the Power button to switch off the machine, then press the Power button again to switch on and restart.

► **If the Start-Up screen stays on:**

- Press the Power button to switch off the machine, then press the Power button again to switch on and restart.
- Contact Help Desk.

► **If the SmartView Xtend resets:**

- Check all the connections and make sure that they are plugged in properly.
- Contact Help Desk if the machine keeps resetting.

● **Returning for Service**

If you have a problem with your SmartView Xtend, call your service agent. Have the machine nearby.

If the problem cannot be fixed over the phone, send your SmartView Xtend to the service agent with a letter describing the problem. If possible, you should use the original SmartView Xtend packaging to avoid damage in transit.

When packing up the machine, make sure that:

- You lock the reading table in its central position by tightening the Margin Stops.
- Your insurance policy covers the SmartView Xtend.

Any SmartView Xtend part found to be defective through normal use within the warranty period will be repaired or replaced free of charge.

You are responsible for freight and insurance charges when you return the equipment to the service agent. However, if your equipment is repaired under warranty, the service agent will cover the return freight charges.

● Specifications

Product

Height: 10"/260 mm

Width: 14"/355 mm

Depth: 17"/430 mm

Weight: 26 lbs/11.60 kg

Operating temperature range: 10 to 40°C

Magnification range

- 15" display: 3x to 37x
- 17" display: 3x to 42x
- 19" display: 4x to 47x
- 21" display: 4x to 52x

Magnification depends on make and model of monitor used.

Illumination

Average lifetime: greater than 10,000 hours

Power Supply

Input voltage: 90-264 V AC

Frequency: 47-63 Hz

Input current: 2.0 A

Power usage: 55 W

PC Compatibility (available with Module 3)

Connects to desktop, laptop, or notebook computer with VGA display. Optional Displays PC and camera images alternately.

Video modes:

- 640x480 @ 60, 72, 75, 85 Hz
- 720x400 @ 70 Hz
- 800x600 @ 60, 72, 75, 85 Hz
- 1024x768 @ 60, 70, 75, 85 Hz
- 1152x864 @ 60, 75 Hz
- 1280x960 @ 60 Hz
- 1280x1024 @ 60 Hz

Optional Accessories

- Smart Keypad upgrade kit
- PC Access upgrade kit
- Simple Remote
- Arm for LCD

Standards and Approvals

See the Declaration of Conformity at the beginning of this manual.

● **Index**

B

Beeps, setting • 18

Brightness button • 10, 13

Brightness, adjusting the • 14

Buttons - Control Panel • 7

C

Calculator • 32, 33

Camera, external • 39, 41

Caring for your Smartview Xtend • 43

Colours

adjusting the image colour • 15

normal/enhanced colour • 16

preferred false colours • 9

setting your false colours • 8

Connections • 2

Contrast, adjusting the • 15

Control Panel • 2, 7

D

Display frequency, adjusting • 16

E

External camera • 39, 41

F

False colours

preferred false colours • 9

setting your false colours • 8

Footswitch • 41

Frequency, adjusting the display • 16

Front panel • See Control Panel

English

G

Getting started • 2, 7

Glare, reducing the • 9, 11

H

Handy hints

- reading • 27

- reducing glare • 9, 11

- setting brightness for text or photos • 11

- signing documents • 23

- using Preset for reading • 24

- viewing cluttered pictures • 26

Help, getting • 45, 51

K

Keypad • 29, 39

- Calculator • 32, 33

- Page key • 23, 34

- PC key • 34, 39

- Size key • 33

- Time and Date • 31, 32

L

Language, setting • 17

Lines and Blinds

- Arrow keys • 35

- Blinds key • 34

- Lines and Blinds button • 24

- Lines key • 34

M

Magnification, adjusting

- See Zoom knob. See Preset key. See Preset button

Main menu • 12

- adjusting the brightness • 14

- adjusting the contrast • 15

- adjusting the image colour • 15

- changing the display frequency • 16

- exiting the menu • 14

Preferences menu • 19,
35, 40
selecting your language
• 17
selecting your support
region • 17
switching audible beeps
on/off • 18
switching between
normal and enhanced
colour • 16
viewing system
information • 18
Maintenance • 43, 51
Margin Stops
location • 2
setting • 11
Menu • 13. See Main
menu
Modules - description • 5
O
On/Off switch • 8

P
Page Locator button • 23
PC key • 34, 39
Platter • See Reading
table
Plus and Minus buttons •
25
Power button • 8
Preferences menu • 19,
35, 40
Preset button • 24
Preset key • 24, 34
R
Reading table • 2
setting Margin Stops •
11
Red light or dot • See
Page Locator button
Region, setting • 17
Remote Control • 21
Lines and Blinds button
• 24

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**Copyright © 2006 der HumanWare Limited,
**Copyright © 2006 van HumanWare Limited,